| CT Injectors |
|--|
| MR Injectors |
| Angio, Vascular and Interventional Injectors 155 |













Look for this symbol throughout the 2013/2014 NEW Diagnostic Imaging Accessories & Supplies catalog for new items this year.

CT Injectors

GE-Nemoto Dual-Shot NCOM Injector

Improve workflow – and save time – with the injector that communicates with your GE CT scanner.

The GE-Nemoto Dual-Shot NCOM Injector is a dual syringe injector designed with GE customer feedback to provide intuitive operation from routine protocols to the most advanced CTA studies. Sold exclusively by GE, this CT injector is integrated with select GE CT scanners for seamless communication. As a result, the potential for technical errors is reduced.

Simply put, the GE-Nemoto Dual-Shot NCOM Injector can help optimize workflow while providing high performance and imaging quality to patients and clinicians alike.

Integrated Communication

- The GE-Nemoto NCOM Injector and GE CT Xtream Injector feature support the CAN CiA 425 communication protocol
- Enables consistent, simultaneous start of injection and scan acquisition protocols

Ease of Use

- Intuitive anatomical protocol selection resembles the GE CT scanner protocol
- User can easily select, modify or create protocols
- Easy syringe loading and unloading no injector command required

Innovation

- Plunger-style syringe allows user to fill syringes outside of the injector, helping to improve patient set-up and throughput
- Quick Fill and Quick Purge functionality is standard
- Cardiac and other advanced protocols, including dilution (simultaneous saline and contrast), multi-phase and saline flush
- Integrated syringe heater enhances patient comfort and lowers injection pressures
- Compact flash memory card allows for quick and easy saving, retrieval, sharing and updating of protocols



Safety Features

- Needle tip pressure and injection parameters displayed in real-time with easy-to-read graphics
- If pressure limit is reached, flow rate decreases automatically and an alert is displayed
- Needle placement test checks patency of the vein to enhance injection safety
- Injector self-tests at power-up and during the entire injection protocol

Syringes and Supplies Available Through GE

- Prepackaged, budget-friendly syringe kits with supplies needed for an exam
- Single and dual syringe kits with Y-tubing and filling spikes or J-tubes

Service

 Backed by GE Healthcare's exceptional service team, providing installation, service and applications training as part of CT scanner training

Upgradeable

• Innovative flash memory card design allows for software changes and upgrades

"One Stop-Shop"

 GE Healthcare provides the complete solution from CT scanner, to injectors, to syringes, and service



| Cat# | Description | List Price | Warranty |
|-----------|--|-------------|----------|
| Injectors | | | |
| E8004M | GE-Nemoto Dual Shot NCOM (Class 4) CT Injector - Pedestal Mount | \$40,600.00 | Е |
| E8004MA | GE-Nemoto Dual Shot NCOM (Class 4) CT Injector - Ceiling Mount for Short Post | \$44,150.00 | Е |
| E8004MB | GE-Nemoto Dual Shot NCOM (Class 4) CT Injector - Ceiling Mount for Medium Post | \$44,150.00 | Е |
| E8004MC | GE-Nemoto Dual Shot NCOM (Class 4) CT Injector - Ceiling Mount for Long Post | \$44,150.00 | Е |
| E8004NG | Nemoto Mounting Plate for Ceiling-Mount Injector | \$650.00 | В |
| | | | |
| GE-Nemoto | Dual-Shot NCOM Syringes and Supplies | | |
| E8004NB | GE-Nemoto Dual Syringe Pack: 200 ml and 100 ml Syringe, Spike and Y-Tube | \$325.00 | Е |
| E8004NC | GE-Nemoto Dual Syringe Pack: 200 ml and 100 ml Syringe w/J-Tube and Y-Tube (20/box) | \$310.00 | Е |
| E8004ND | GE-Nemoto Single Syringe Pack: 200 ml Syringe w/Spike and 60 in. Coiled Tubing (50/box) | \$480.00 | Е |
| E8004NE | GE-Nemoto Single Syringe Pack: 200 ml Syringe w/J-Tube and 60 in. Coiled Tubing (50/box) | \$465.00 | Е |
| E8004RA | GE-Nemoto 200 ml Syringe w/J-Tube (50/Case) | \$375.00 | Е |
| E8004RB | GE-Nemoto 100 ml Syringe w/J-Tube (50/Case) | \$340.00 | Е |
| E8004RC | GE-Nemoto J-Tube only (100/Case) | \$135.00 | Е |
| E8004RD | GE-Nemoto 60 in. Straight Line (100/Case) | \$162.00 | Е |
| E8004RE | GE-Nemoto 60 in. Y-Tube, 1 Check Valve, Straight (100/Case) | \$415.00 | Е |
| E8004RF | GE-Nemoto 200 ml Syringe only (50/Case) | \$370.00 | Е |
| E8004RG | GE-Nemoto 100 ml Syringe only (50/Case) | \$330.00 | Е |
| E8004RH | GE-Nemoto Fill Set (100/Case) | \$470.00 | Е |
| E8004RJ | GE-Nemoto 60 in. Coiled Line (100/Case) | \$225.00 | Е |
| E8004RK | GE-Nemoto 60 in. Coiled Line Y-Tube, 1 Check Valve (100/Case) | \$430.00 | Е |
| E8004RL | GE-Nemoto 60 in. Y-Tube, 2 Check Valves, Straight (100/Case) | \$475.00 | Е |
| E8004RM | GE-Nemoto 200/200 ml Syringe w/ spike (20/box) | \$395.00 | Е |
| E8004RN | GE-Nemoto 200/200 ml Syringe w/ fill tube (20/box) | \$383.00 | Е |
| E8004RP | GE-Nemoto Fill Set, 7 in., with spike, cap and clamp (100/Case) | \$470.00 | Е |

GE-Nemoto CT Dual-Shot Injection System

Brought to you exclusively by GE Healthcare

The GE-Nemoto Dual-Shot Injector rounds out the GE Accessories family of products for your complete CT suite. Sold exclusively by GE, this CT injector helps optimize contrast procedures and contrast usage as well as streamline clinical workflow. GE-Nemoto ... It's That Simple!

Ease of Use

- Intuitive anatomical protocol selection complements the GE CT scanner protocol set-up
- Quick and easy retrieval and programming with its intuitive touch screen
- Easy syringe loading and unloading
- Winner of Good Design Award

Innovation

- Innovative plunger style syringe allows user to fill syringes outside the injector, improving patient set-up
- Quick Fill and Quick Purge functionality is standard
- Advanced protocols with dilution (simultaneous saline and contrast), multi-phases and saline flush
- Compact flash card for protocol saving, quick retrieval, sharing and easy updates

Syringes and Supplies

- Prepackaged budget-friendly syringe kits with supplies needed for an exam
- Single and dual syringe kits for contrast only and saline exams with filling option J-tube or spike and Y-tubing
- GE will work with you to select the appropriate products, create a blanket purchase order and schedule regular deliveries of any configuration and quantity of products that you require





Safety Features

- Real-time pressure monitoring graphically monitors needle tip pressure and injection parameters, as well as a needle placement test to help administer injections
- Injector self-tests at power-up and during the entire injection protocol

Service

 Backed by GE Healthcare's exceptional service team, providing installation, service and applications training as part of CT scanner training

Upgradeable

• Flash memory card design allows for software changes and upgrades

"One-Stop-Shop"

• GE Healthcare provides the complete solution from CT scanner, to injectors, to syringes, and service

| Cat# | Description | List Price | Warranty |
|--------------|--|-------------|----------|
| Injectors | | | |
| E8004N | GE-Nemoto Dual-Shot CT Injector Pedestal-Mount | \$34,000.00 | Е |
| E8004NA | GE-Nemoto Dual-Shot CT Injector Ceiling-Mount (Short Post – ceiling height: 8 ft. 2 in. – 8 ft. 10 in.) REQUIRES E8004NG Mounting Plate be added to order | \$37,500.00 | E |
| E8004NJ | GE-Nemoto Dual-Shot CT Injector Ceiling-Mount (Medium Post – ceiling height: 9 ft. 0 in. – 10 ft. 6 in.) REQUIRES E8004NG Mounting Plate be added to order | \$37,500.00 | Е |
| Injector Acc | essories | | |
| E8004NG | Nemoto Mounting Plate for Ceiling-Mount Injector | \$650.00 | В |



GE-Nemoto Dual-Shot NCOM Syringes and Supplies


| | , , , | | |
|---------|--|----------|---|
| E8004NB | GE-Nemoto Dual Syringe Pack: 200 ml and 100 ml Syringe, Spike and Y-Tube | \$325.00 | Е |
| E8004NC | GE-Nemoto Dual Syringe Pack: 200 ml and 100 ml Syringe w/J-Tube and Y-Tube (20/box) | \$310.00 | Е |
| E8004ND | GE-Nemoto Single Syringe Pack: 200 ml Syringe w/Spike and 60 in. Coiled Tubing (50/box) | \$480.00 | Е |
| E8004NE | GE-Nemoto Single Syringe Pack: 200 ml Syringe w/J-Tube and 60 in. Coiled Tubing (50/box) | \$465.00 | Е |
| E8004RA | GE-Nemoto 200 ml Syringe w/J-Tube (50/Case) | \$375.00 | Е |
| E8004RB | GE-Nemoto 100 ml Syringe w/J-Tube (50/Case) | \$340.00 | Е |
| E8004RC | GE-Nemoto J-Tube only (100/Case) | \$135.00 | Е |
| E8004RD | GE-Nemoto 60 in. Straight Line (100/Case) | \$162.00 | Е |
| E8004RE | GE-Nemoto 60 in. Y-Tube, 1 Check Valve, Straight (100/Case) | \$415.00 | Е |
| E8004RF | GE-Nemoto 200 ml Syringe only (50/Case) | \$370.00 | Е |
| E8004RG | GE-Nemoto 100 ml Syringe only (50/Case) | \$330.00 | Е |
| E8004RH | GE-Nemoto Fill Set (100/Case) | \$470.00 | Е |
| E8004RJ | GE-Nemoto 60 in. Coiled Line (100/Case) | \$225.00 | Е |
| E8004RK | GE-Nemoto 60 in. Coiled Line Y-Tube, 1 Check Valve (100/Case) | \$430.00 | Е |
| E8004RL | GE-Nemoto 60 in. Y-Tube, 2 Check Valves, Straight (100/Case) | \$475.00 | Е |
| E8004RM | GE-Nemoto 200/200 ml Syringe w/ spike (20/box) | \$395.00 | Е |
| E8004RN | GE-Nemoto 200/200 ml Syringe w/ fill tube (20/box) | \$383.00 | Е |
| E8004RP | GE-Nemoto Fill Set, 7 in., with spike, cap and clamp (100/Case) | \$470.00 | Е |

GE-Nemoto Body Weight Protocol

Help reduce protocol variation and help improve workflow.

GE-Nemoto injector's alternative programming option provides medical professionals the opportunity to enter known information such as: patient weight, contrast media concentration, desired iodine amount, injection time and volume. The Body Weight Protocol then calculates and displays the equivalent injection protocol in terms of flow rate and total volumes.

- Protocol calculated based upon patient weight and desired iodine amount
- Helps reduce contrast waste
- Helps improve protocol consistency within a department
- Customizable interface to meet the needs of the user
- Helps improve injection performance from patient to patient

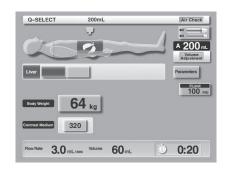

| Cat# | Description | List Price | Warranty |
|---------|--------------------------------|------------|----------|
| E8004NL | GE-Nemoto Body Weight Protocol | \$5,000.00 | Е |

MEDRAD Stellant® D with DualFlow Injection Systems

Stellant D with DualFlow is more than a saline flush after a contrast bolus. You can inject both contrast and saline at the same time. The key is the simultaneous injection capability of the DualFlow option. DualFlow enables variable ratios of plunger motion from contrast and saline syringes simultaneously. With the proper ratio, left and right heart ventricles can be illuminated uniformly for good image quality.

- Stellant D with DualFlow incorporates all the benefits of Stellant D. It also allows users to test vein accesses with saline and prime patient tubing with saline to save contrast
- Benefits for Stellant D with DualFlow:
 - Helps match attenuation levels in both ventricles
- Additional benefits include:
 - Automation features to help maximize throughput: integrated auto load, auto retract, auto prime and auto syringe sensing

- Enhanced syringe design for quick, easy install and detachment
- Remote check for air confirmation button on injector head; remote arming at injector head
- Pressure monitor graph and flow profile preview
- Up to 6 phases, including pause and hold capabilities
- Programmable pressure limit
- Color touch screen
- Available in either ceiling counterpoise or pedestal-mount configurations



All Ceiling-Mount MEDRAD Injectors require catalog #E8007PJ.

| Cat# | Description | List Price | Warranty |
|----------------|--|-------------|----------|
| Injectors | | | |
| E8007NH | Pedestal-Mount Injection System with DualFlow | \$42,000.00 | Е |
| E8007NJ | Ceiling-Mount (Long Post – ceiling height: 10 ft. 6 in 12 ft.) Injection System w/DualFlow REQUIRES E8007PJ Mounting Plate be added to the order | \$47,500.00 | Е |
| E8007PP | Ceiling-Mount (Medium Post – ceiling height: 9 ft. 0 in. – 10 ft. 6 in.) Injection System w/DualFlow REQUIRES E8007PJ Mounting Plate be added to the order | \$47,500.00 | Е |
| E8007NG | Ceiling-Mount (Short Post – ceiling height: 8 ft. 2 in. – 8 ft. 10 in.) Injection System w/DualFlow REQUIRES E8007PJ Mounting Plate be added to the order | \$47,500.00 | Е |
| E8007NW | Stellant D w/DualFlow Mobile CT Injection System | \$51,500.00 | Е |
| Injector Acces | ssories | | |
| E8007PJ | Mounting Plate for Overhead Counterpoise System (OCSIII) | \$650.00 | В |
| E8007PG | Stellant OCS Installation by MEDRAD | \$3,500.00 | Е |
| E8007PH | Stellant Pedestal Injector Installation by MEDRAD | \$2,600.00 | Е |
| Syringes and | Supplies | | |
| E8006RB | Stellant D & D w/DualFlow Dual Syringe Quad-Pak w/Spike (20/box) (SDS-CTP-SPK) | \$535.00 | Н |
| E8006RC | Stellant D & D w/DualFlow Dual Syringe Quad-Pak w/Quick Fill Tube (20/box) (SDS-CTP-QFT) | \$485.00 | Н |
| E8006RD | Stellant Sx* D & D w/DualFlow Swabbable-Valve Transfer Set (50/box) (SVTS) | \$475.00 | Н |



MEDRAD Stellant D Injection Systems

- Stellant D features saline flush capability. It also allows users to test vein accesses with saline; prime patient tubing with saline to save contrast
- Both models include:
- Automation features to help maximize throughput: integrated auto load, auto retract, auto prime and auto syringe sensing
- Enhanced syringe design for quick, easy install and detachment
- Remote check for air confirmation button on injector head; remote arming at injector head
- Pressure monitor graph and flow profile preview
- Up to 6 phases including pause and hold capabilities
- Programmable pressure limit
- Color touch screen
- Available in either ceiling counterpoise or pedestal-mount configurations



All Ceiling-Mount MEDRAD Injectors require catalog #E8007PJ.

| Cat# | Description | List Price | Warranty |
|---------------|--|-------------|----------|
| Injectors | | | |
| E8007NC | Pedestal-Mount Injection System | \$36,750.00 | Е |
| E8007NF | Ceiling-Mount (Long Post – ceiling height: 10 ft. 6 in 12 ft.) Injection System REQUIRES E8007PJ Mounting Plate be added to the order | \$42,000.00 | E |
| E8007PN | Ceiling-Mount (Medium Post – ceiling height: 9 ft. 0 in. – 10 ft. 6 in.) Injection System w/DualFlow REQUIRES E8007PJ Mounting Plate be added to the order | \$42,000.00 | E |
| E8007ND | Ceiling-Mount (Short Post – ceiling height: 8 ft. 2 in. – 8 ft. 10 in.) Injection System REQUIRES E8007PJ Mounting Plate be added to the order | \$42,000.00 | E |
| E8007NT | Stellant D Mobile CT Injection System | \$45,220.00 | Е |
| Injector Acce | essories | | |
| E8007PJ | Mounting Plate for Overhead Counterpoise System (OCSIII) | \$650.00 | В |
| E8007PG | Stellant OCS Installation by MEDRAD | \$3,500.00 | Е |
| E8007PH | Stellant Pedestal Injector Installation by MEDRAD | \$2,600.00 | Е |
| E8007NP | Stellant D DualFlow Upgrade Kit | \$7,100.00 | Е |
| Syringes and | Supplies | | |
| E8006RB | Stellant D & D w/DualFlow Dual Syringe Quad-Pak w/Spike (20/box) (SDS-CTP-SPK) | \$535.00 | Н |
| E8006RC | Stellant D & D w/DualFlow Dual Syringe Quad-Pak w/Quick Fill Tube (20/box) (SDS-CTP-QFT) | \$485.00 | Н |
| E8006RD | Stellant Sx D & D w/DualFlow Swabbable-Valve Transfer Set (50/box) (SVTS) | \$475.00 | Н |

MEDRAD Stellant Injection System with Integrated CT Communication

GE's Enhanced Xtream Injector option can now be coupled with MEDRAD's Stellant CT Injector to provide clinicians with a simplified and integrated imaging solution. The result is simplified CT scanning that enables improved bolus timing and helps optimize workflow.

Intelligent, Intuitive and Innovative

MEDRAD's Stellant D is the smart choice for dependable performance exactly when you need it – during critical CT contrast exams, where accuracy is vital and real-time injection graphics are indispensable. The Stellant D is consistently reliable, easy to use, and offers enhanced safety features, from the dual-syringe capability, automatation features, to the Enhanced Xtream injector interface.

Integrated Communication Advantages

The MEDRAD Stellant CT Injection System is validated for use with GE's Enhanced Xtream Injector option on selected scanners – enabling CAN Class 4 functionality for seamless communication. The resulting injector and CT scanner integration benefits include:

- Reduced overall programming time
- Improved scanner and injector protocol matching through programming of the injector from the scanner console
- Better control of contrast enhancement, since the contrast injection and CT scan start time are synchronized
- Efficient documentation of the injection results in PACS with automatic transfer of contrast injection results



| Cat# | Description | List Price | Warranty |
|--------------|--|-------------|----------|
| E8037NH | MEDRAD Stellant D Dualflow Pedestal Mount ISI Ready with Integrated CT Communication | \$52,625.00 | E |
| E8037NG | MEDRAD Stellant D Dualflow Ceiling Mount for Short Post ISI Ready with Integrated CT Communication | \$58,125.00 | Е |
| E8037PP | MEDRAD Stellant D Dualflow Ceiling Mount for Medium Post ISI Ready with Integrated CT Communication | \$58,125.00 | Е |
| E8037NJ | MEDRAD Stellant D Dualflow Ceiling Mount for Long Post ISI Ready with Integrated CT Communication | \$58,125.00 | Е |
| For Upgrades | | | |
| E8007CB | MEDRAD Stellant Integrated Injector ISI 900 (Upgrades for E8007NG, E8007NH, E8007NJ, E8007PP Injectors Only) | \$10,625.00 | E |

MEDRAD XDS™ Extravasation Detector Technology for the Stellant CT Injection System

- Directly senses contrast or saline pooling under a patient's skin during a CT procedure
- XDS uses radiofrequency (RF) waves to reliably detect an extravasation and stop the injection
- The patented XDS technology is designed to prevent moderate to severe extravasations
- XDS Patient Interface Disposables are designed to allow easy and consistent positioning of the XDS sensors



| Cat# | Description | List Price | Warranty |
|---------|--|-------------|----------|
| E8007NK | MEDRAD XDS Extravasation Detector Technology | \$16,000.00 | Н |



MEDRAD Stellant Sx Injection Systems

- Stellant Sx is a single, contrast-only CT injection system designed for fast-paced CT sites that require maximum patient throughput
- Models include:
 - Automation features to help maximize throughput: integrated auto load, auto retract, auto prime and auto syringe sensing
 - Enhanced syringe design for quick, easy install and detachment
 - Remote check for air confirmation button on injector head; remote arming at injector head
 - Pressure monitor graph and flow profile preview
 - Up to 6 phases, including pause and hold capabilities
 - Programmable pressure limit
 - Color touch screen
 - Available in either ceiling counterpoise or pedestal-mount configurations




All Ceiling-Mount MEDRAD Injectors require catalog #E8007PJ.

| Cat# | Description | List Price | Warranty |
|---------------|--|-------------|----------|
| Injectors | | | |
| E8007NA | Pedestal-Mount Injection System | \$28,000.00 | Е |
| E8007NE | Ceiling-Mount (Long Post – ceiling height: 10 ft. 6 in. – 12 ft.) REQUIRES E8007PJ Mounting Plate be added to the order | \$33,000.00 | E |
| E8007PK | Ceiling-Mount (Medium Post – ceiling height: 9 ft. 0 in. – 10 ft. 6 in.) Injection System w/DualFlow REQUIRES E8007PJ Mounting Plate be added to the order | \$33,000.00 | E |
| E8007NB | Ceiling-Mount (Short Post – ceiling height: 8 ft. 2 in. – 8 ft. 10 in.) Injection System REQUIRES E8007PJ Mounting Plate be added to the order | \$33,000.00 | Е |
| E8007NS | Stellant Sx Mobile CT Injection System | \$36,000.00 | Е |
| Injector Acce | essories | | |
| E8007PJ | Mounting Plate for Overhead Counterpoise System (OCSII) | \$650.00 | В |
| E8007NR | Stellant Sx to D Dual Head Upgrade Kit | \$15,500.00 | Е |
| E8007PG | Stellant OCS Installation by MEDRAD | \$3,500.00 | Е |
| E8007PH | Stellant Pedestal Injector Installation by MEDRAD | \$2,600.00 | Е |
| Syringes and | Supplies | | |
| E8006RE | Stellant Sx Single Syringe Tri-Pak w/Spike (50/box) (SSS-CTD-SPK) | \$885.00 | Н |
| E8006RF | Stellant Sx Single Syringe Tri-Pak w/Quick Fill Tube (50/box) (SSS-CTP-QFT) | \$825.00 | Н |
| E8006RD | Stellant Sx D & D w/DualFlow Swabbable-Valve Transfer Set (50/box) (SVTS) | \$475.00 | Н |

MEDRAD P3T® for Stellant D CT Injectors

- Enables generation of individualized contrast injection protocols, customized to the patient
- Utilizes the power of DualFlow technology (simultaneous injection of contrast and saline) to obtain functional cardiac data
- Helps provide consistency amongst clinicians
- Used with CT angiography of cardiac structures, coronary arteries, chambers of the heart, thoracic and abdominal aorta, and pulmonary vasculature
- Available on the MEDRAD Stellant D CT Injection System

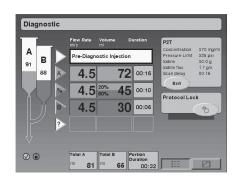

| Cat# | Description | List Price | Warranty |
|--------------------|--|------------|----------|
| NEW E8007PT | MEDRAD P3T Cardiac Option | \$5,300.00 | Н |
| NEW E8007PS | MEDRAD P3T Abdomen Option | \$5,500.00 | Н |
| NEW E8007BA | MEDRAD P3T PA Pulmonary Angiography Option | \$5,500.00 | Н |

MEDRAD Connect PACS

Connect.PACS is a radiology informatics application that integrates each patient's CT contrast injection record with the associated clinical image set stored in PACS.

- Communicates to radiologists contrast injection protocol details associated with the CT images stored on PACS
- Enables radiologists to improve protocols by retrospectively analyzing contrast injection protocols along with the associated image set
- Automatically records consistent data throughout all shifts, and eliminates the need to manually enter or scan contrast information into PACS
- Facilitates and simplifies the processes associated with Joint Commission compliance
- Reduces discrepancies in reconciling dictated patient records with contrast-enhanced CT billing and record keeping
- Facilitates robust reporting and analysts to identify new improvement opportunities such as reducing contrast waste, increasing productivity and reducing atypical injection events
- AutoCapture instantly transfers contrast injection information from the Stellant CT Injector

- AutoMatch integrates injection information with study information including patient ID, name and accession number available from the scanner
- ArchiveTransfer automatically pushes injection information to the PACS archive
- Access enables visibility into aggregated contrast injection information from any PC connected to your hospital network
- Robust Reporting capabilities provide visibility into contrast injection and management metrics as well as atypical injection events
- Data Export enables data to be exported into spreadsheets for more detailed analysis and reporting
- Print options include the ability to print to any standard printer on the network, or to a label printer
- Security and compliance safeguards are provided through features including authentication, access control and audit control

| | Cat# | Description | List Price |
|-----|---------|---------------------|-------------|
| NEW | E8007BB | MEDRAD Connect.PACS | \$36,000.00 |



ACIST® Empower CTA Double Barrel Injector

The ACIST CTA injector provides the ability to inject a saline chase and also has a patented extravasation detection accessory (EDA). This may help to protect patients from extravasations (E8005JD and E8005JE).

Design Enhanced Productivity

- IRiSCT Contrast Management tool. Automatically captures and consolidates essential contrast and syringe usage data
- Full rotation of injector head, either to the right or left a full 270,° allows use from either side of the gantry or table
- The injector controls swivel and remains vertical at all times. The LED display is bright with high contrast and easy-to-see screen even at sharp angles
- Well-labeled buttons and matching lights color-coded syringe identification system
- Arming at the injector enhances efficiency and throughput by allowing operator to arm injector in CT room. Great feature for small hospitals, imaging centers or facilities that operate CT rooms with one individual

Safety Features

- Patented EDA helps detect estravasation events
- Tilt sensor/lockout helps reduce possible injection of air embolus
- Dual pressure displays at injector head and remote for easy monitoring
- Voice prompts announce status of injection details for operator



| Cat# | Description | List Price | Warranty |
|---------|-------------------------------|-------------|----------|
| E8005JF | Pedestal-Mount Injector | \$34,750.00 | Е |
| E8005JG | Ceiling-Mount Injector | \$38,750.00 | Е |
| E8005JD | Pedestal-Mount Injector w/EDA | \$36,750.00 | Е |
| E8005JE | Ceiling-Mount Injector w/EDA | \$40,750.00 | Е |

ACIST Empower CT Injector

- Extravasation detection: ACIST Extravasation Detection Accessory (EDA) technology can help detect extravasations
- Pendant Switch: Puts convenient "on/off" procedure control in your hand
- Injector Knob: Provides an easy means of manually advancing the syringe to eliminate air from the line and for checking patient patency; rotates to offer a highly visible indication that injection is in progress
- Tilt Sensor/Lockout: Injector head tilts in either direction for easy positioning prior to injection; helps minimize air embolism risk
- Syringe Fill Options: Choose auto or manual syringe filling or preload syringes with your choice of contrast
- Injector Controller: Syringe filling, removal and monitoring even system arming and flow rate manipulation can be controlled at the injector head



| Cat# | Description | List Price | Warranty |
|---------|--|-------------|----------|
| E8005J | ACIST Empower CT Injector – Pedestal-Mount | \$20,000.00 | E |
| E8005JB | ACIST Empower CT Injector w/EDA – Pedestal-Mount | \$21,500.00 | Е |

E-Z-EM™ Protocol™ Colon Insufflator for Virtual Colonoscopy

E-Z-EM has an Automated Carbon Dioxide Insufflation System for Virtual Colonoscopy.

- Patient Comfort Carbon dioxide is readily absorbed into the bowel and may improve
 patient comfort after the procedure, compared with room air insufflation. Helps
 ensure patient safety through redundant pressure relief valves and automatic flow
 stop features, which protects against over-insufflation. A choice of tips for rectal
 administration allows the practitioner to select the most comfortable tip for the
 patient
- Consistent Insufflation Automated insufflation to a constant pressure helps take the
 guesswork out of determining the endpoint. Regulated pressure gently distends the
 colon. Volume display allows the monitoring and recording of the volume of carbon
 dioxide used, unlike room air insufflation, where it is impossible to judge the amount
 of air introduced
- Productivity Automated insufflation may decrease the demands on staff time
- E8005JL package includes the PROTOCOL Colon Insufflator Accessory Cart



| Cat# | Description | List Price | Warranty |
|---------|---|-------------|----------|
| E8005JL | E-Z-EM PROTOCOL Colon Insufflator Package | \$11,349.00 | E |



MR Injectors

MEDRAD Spectris® Solaris EP MR Injector with Integrated Continuous Battery Charger

- Injector:
- Designed for use in all MR scanner field strengths up to and including 3T
- Longer battery life with increased capacity
- Integrated battery charger
- Independent KVO (Keep Vein Open) allows patient to be prepared before beginning the scan and to move from noncontrasted to contrasted scan
- Larger 115 ml saline syringe for longer KVO or multiple flushes
- Optimized touch screen with fewer keystrokes
- Six user-programmable phases for added flexibility
- Fiber optic cable for direct and reliable communication
- Mobile Injector:
- Same benefits and features as the Solaris EP, validated for the MR mobile imaging environment
- Solaris EP Accessories:
- Integrated Continuous Battery Charger (ICBC)
- Syringe loader
- Handswitch



- MRI Disposable Kit for MEDRAD Solaris EP Injector consists of:
 - (2) FLS syringes with a fluid delivery device
 - (1) 115 ml Qwik-Fit Syringe* for saline
 - (1) 65 ml Qwik-Fit Syringe for contrast
 - (1) large spike, (1) small spike
 - (1) 96 in. low pressure connector tube with T-connector and check valve
 - Packaged 50 disposable kits per case

| Cat# | Description | List Price | Warranty |
|-----------------|--|-------------|----------|
| Injectors | | | |
| E8804SB | MEDRAD Spectris Solaris EP Injection System | \$49,375.00 | Е |
| E8804SN | MEDRAD Spectris Solaris EP Injection System for Mobile | \$57,000.00 | Е |
| Injector Access | sories | | |
| E8804SZ | MR Solaris EP ICBC | \$2,700.00 | Н |
| E8804SR | MEDRAD Solaris EP Installation | \$2,900.00 | N/A |
| Syringes and S | Supplies | | |
| E8804SA | MRI Disposable Kit for MEDRAD Solaris (50/case) | \$1,400.00 | Н |

SmartSet for MRI or MRA Contrast Injection TopSpins Handheld MR Injection System SmartSet (Side Port Availability)

- Enhances the reliability of bolus timing
- Helps minimize the bolus gap during the switch of gadolinium to saline flush
- Tubing designed specifically for the dynamic gadolinium injections of MRA or MRI
- Swivel Lock to help prevent IV separation from angio catheter
- Two inputs for simultaneous attachment of contrast and flush syringes
- Check valves automatically switch between contrast and flush
- Kink resistant, adjustable length, transparent, latex-free

TopSpins MR Adapter

- Connects to standard IV extension tubing of any length
- Luer Lock to help prevent separation from IV tubing



| Cat# | Description | List Price | Warranty |
|---------|---|------------|----------|
| E8804LA | TopSpins SmartSet | \$900.00 | Н |
| E8804LB | TopSpins MR Adapter | \$500.00 | Н |
| E8804LC | TopSpins SmartSet w/Side Port (15/case) | \$450.00 | Н |
| E8804LD | TopSpins SmartSet w/Side Port (30/case) | \$825.00 | Н |



Angio, Vascular and Interventional Injectors

MEDRAD Mark 7 Arterion Injector

The Mark 7 Arterion™ is light, maneuverable and easy to use. Less time positionir and setting up the Arterion means more time with the patient. The clearly visible and intuitive user interface guides you through proper set-up, and highlights the information you need to perform injections confidently.

- Ergonomic handle for easier maneuverability
- Front load syringe for simple insertion and clean removal
- Syringe provides a clear view of the contrast
- Light injector head with a handle to make it easier to position for injection
- Smooth arc design of pedestal for extended reach
- Small footprint which increases mobility around a busy lab
- Bright and colorful intuitive user interface is designed to highlight the information you need
- Highlighted armed state to know when the system is ready to inject
- History and protocol screens to easily access the amount of contrast delivered to the patient and store and recall protocols



| | Cat# | Description | List Price | Warranty |
|-----|---------|--|-------------|----------|
| NEW | E7018HA | Mark 7 Arterion injector on integrated pedestal with installation and warranty (GE Innova package) | \$45,000.00 | E |
| NEW | E7018HC | Mark 7 Arterion injector on integrated pedestal with installation and warranty (GE OEC package) | \$45,000.00 | Е |
| NEW | E7018HD | Mark 7 Arterion injector on table mount with installation and warranty (GE OEC package) | \$47,000.00 | Е |

ACIST® CVi® Contrast Delivery System

An advanced variable rate contrast injection system developed for use in cardiac and vascular angiography procedures, including coronary, ventricle, aorta and peripheral vasculatures. The CVi system is designed to replace the lab's hand manifold and power injector which means efficiency and ease of operation for the physician and staff. The responsive variable rate delivery of contrast provides quality imaging with less contrast used and delivered to the patient.

Use of the system's AngioTouch® hand controller allows real-time, responsive control of contrast delivery – from a puff to a full injection. It also helps minimize radiation dose to the staff because the hand controller allows the operator to step back from the radiation source.



Key features include:

- AngioTouch hand controller facilitates real-time, responsive control of contrast, contributing to reduced contrast dosage while providing quality images
- Air column detection sensor stops injection if air is detected, alerting operator of potential safety concern
- Automatic syringe-refill and air purge reduce operator steps, providing ease of operation and efficiency
- Integrated pressure transducer provides continuous hemodynamic monitoring, automatically sending real-time readings to lab monitors
- Touch screen interface provides easy set up of injection parameter limits, tracking contrast delivered while providing monitoring information and system alerts

| Cat# | Description | List Price | Warranty |
|---------|--|-------------|----------|
| E7016NC | ACIST CVi Contrast Delivery System - Pedestal or Table Mounted | \$36,250.00 | Е |
| E7016NE | ACIST CVi Contrast Delivery System - Pedestal Mounted | \$35,600.00 | Е |
| E7016NH | ACIST CVi Contrast Delivery System - Table Mounted | \$35,600.00 | Е |

MEDRAD Avanta™ Injection Systems

MEDRAD Avanta is an advanced fluid management system that is flexible enough for multiple cardiac imaging applications; coronary, ventriculograms and aortograms, and more. This system is intended to inject contrast in variable mode at low flow rates and low pressures for coronary and small vessels, fixed mode at high flow rates and high-pressure injections for larger vessels and peripherals, and fixed mode for saline flushing. Of course, in MEDRAD tradition, the MEDRAD Avanta has been designed with safety, ease-of-use and reliability in mind.

- Enhanced Safety:
 - Precise control of coronary injection, especially with small catheters
 - Accurate injection pressure control with user-adjustable pressure limits
- Ease-of-use:
 - Easy-to-read, color touch screen with On-screen Tutorial for simplified set-up
 - Multi-patient disposables (up to 5 patients)

- Precision:
 - Distinct bolus sharpness feature
 - Exact variable and fixed-flow delivery with increased accuracy and responsiveness with hand controller (when compared to hand injection)
 - Automatic transition from injection to hemodynamic state
- Flexibility:
 - Pedestal or table mount with independent mounting of display and injection unit



| Cat# | Description | List Price | Warranty |
|---------|--|-------------|----------|
| E7018MA | Avanta Pedestal-Mount Injection System | \$58,000.00 | E |
| E7018MB | Avanta Table-Mount Injector and 90 ft. TIF | \$61,000.00 | Е |



MEDRAD Mark V ProVis® Injection Systems

- Mechanical stop to minimize risk of over-injection
- Articulating arm extends reach for ease in positioning close to patient
- Height-adjustable arm for ease in loading syringes
- Rotating, tilting control panel is easy to see and read from almost any position in the room
- Multiple syringe configurations (60/150, 150/150 or 200/200 ml) help to match the right syringe size with the patient and procedure – choose either 150 or 200 ml dual syringe capacity
- Floor-, ceiling- or wall-mount options add flexibility
- Segmented control panel is easy to learn and operate
- "One Touch" enable button, and variable speed control, for faster loading and easier delivery



All Ceiling-Mount MEDRAD Injectors require catalog #E8007PJ.

| Cat# | Description | List Price | Warranty |
|---------------|---|-------------|----------|
| Injectors | | | |
| E7018JP | Pedestal-Mount Injection System | \$27,000.00 | N/A |
| E7018JR | Table-Mount Injection System w/Remote Keyboard | \$32,000.00 | N/A |
| E7018JT | Ceiling-Mount Injection System REQUIRES E8007PJ Mounting Plate be added to the order | \$36,000.00 | N/A |
| E9200AD | Ceiling-Mount Injection System and Track REQUIRES E8007PJ Mounting Plate be added to the order | \$40,000.00 | N/A |
| E7018LM | Pedestal Injection System for Mobile Applications | \$38,000.00 | N/A |
| E7018JN | Table-Mount Injector w/Keyboard and Free-Standing Pedestal | \$33,500.00 | N/A |
| E7018JW | Movable Ceiling-Mount Injector | \$35,000.00 | N/A |
| Injector Acce | essories | | |
| E7018JS | Interface Cable for OEC | \$1,500.00 | N/A |
| E7018JZ | Mavig 2.5m Track without Cable Spooler | \$3,200.00 | N/A |
| E8007PJ | Mounting Plate for Overhead Counterpoise System (OCSII) | \$650.00 | В |
| Syringes and | d Supplies | | |
| E8006R | 150 ml FasTurn™ Syringe with Tube | \$475.00 | Е |
| E8006RL | 200 ml Syringe QFT | \$600.00 | Е |